


## Omega Optical Filter Set XF02-2

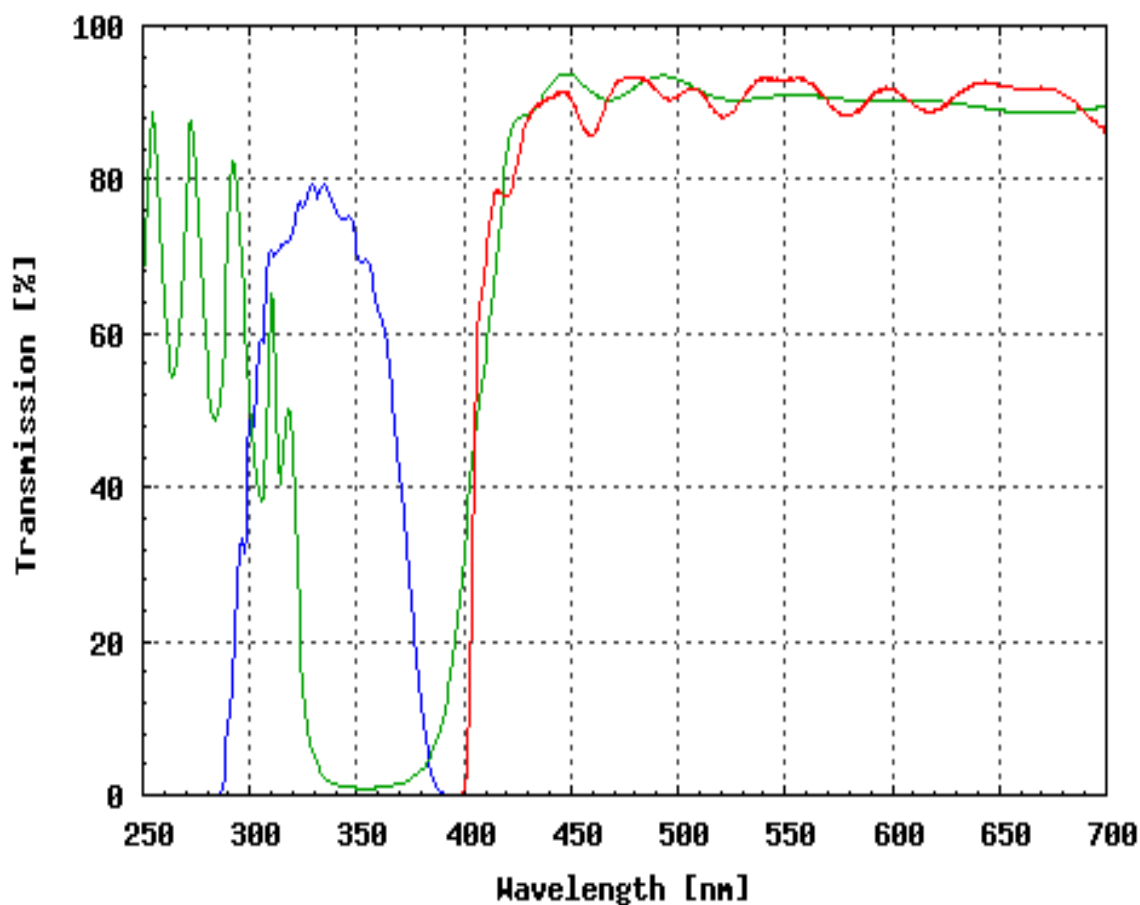

### FILTER SET COMPONENTS

Exciter:	XF1001 (330WB80)
Dichroic:	XF2001 (400DCLP)
Emitter:	XF3097 (400ALP)

### OTHER INFORMATION

Longpass blue set with a 330nm wide-band exciter. Brighter emission with lower signal-to-noise than the XF03. Brightest image and background; however, due to broad excitation filter, may photobleach sample. Colours seen through this filter can be variable depending on exact emission of dye. For use with Xe arc lamp. Can also be used for all Quantum Dot Conjugates, for simultaneous multi-colour viewing with Xenon excitation.