


Omega Optical Filter Set XF03

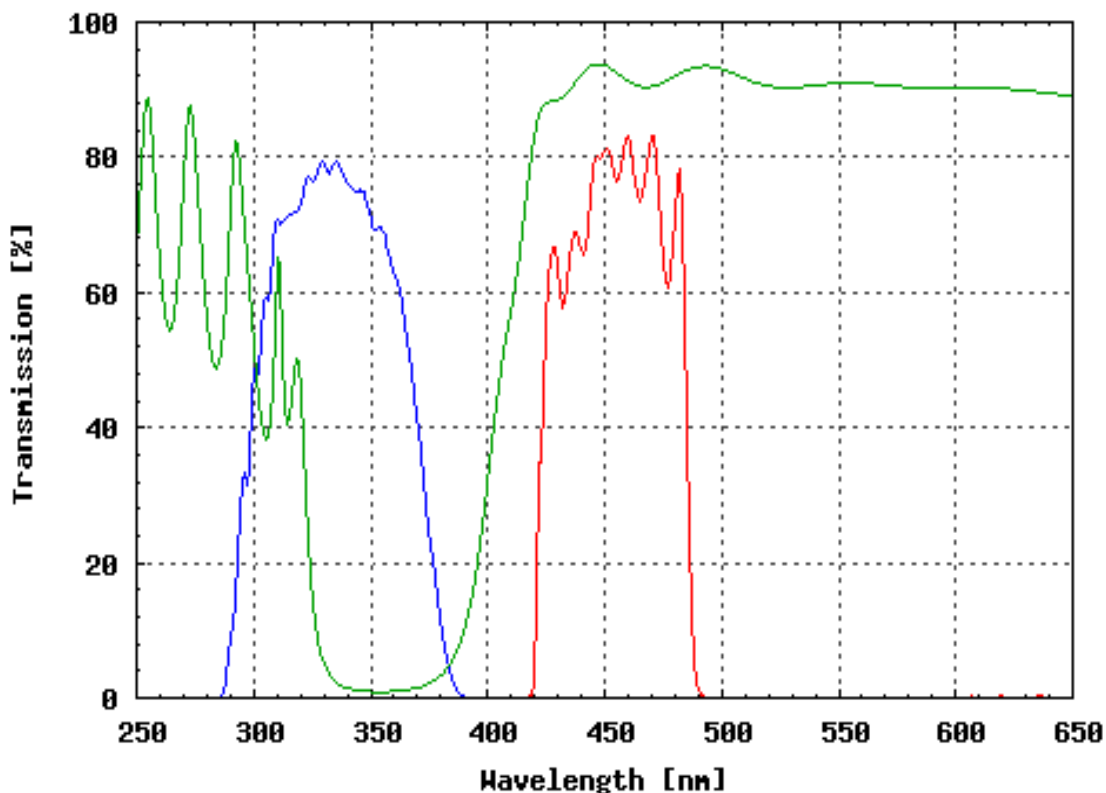

FILTER SET COMPONENTS

Exciter:	XF1001 (330WB80)
Dichroic:	XF2001 (400DCLP)
Emitter:	XF3002 (450DF65)

OTHER INFORMATION

Bandpass blue set with a 330nm wide-band exciter. Highly recommended for multicolour work using a xenon (Xe) light source. This set will also work for a mercury (Hg) light source, but in this case, the set will have a brighter background than our XF06 filter set, which is more optimised for mercury arc lamps. Higher signal-to-noise than the XF02-2.