


## Omega Optical Filter Set XF05-2

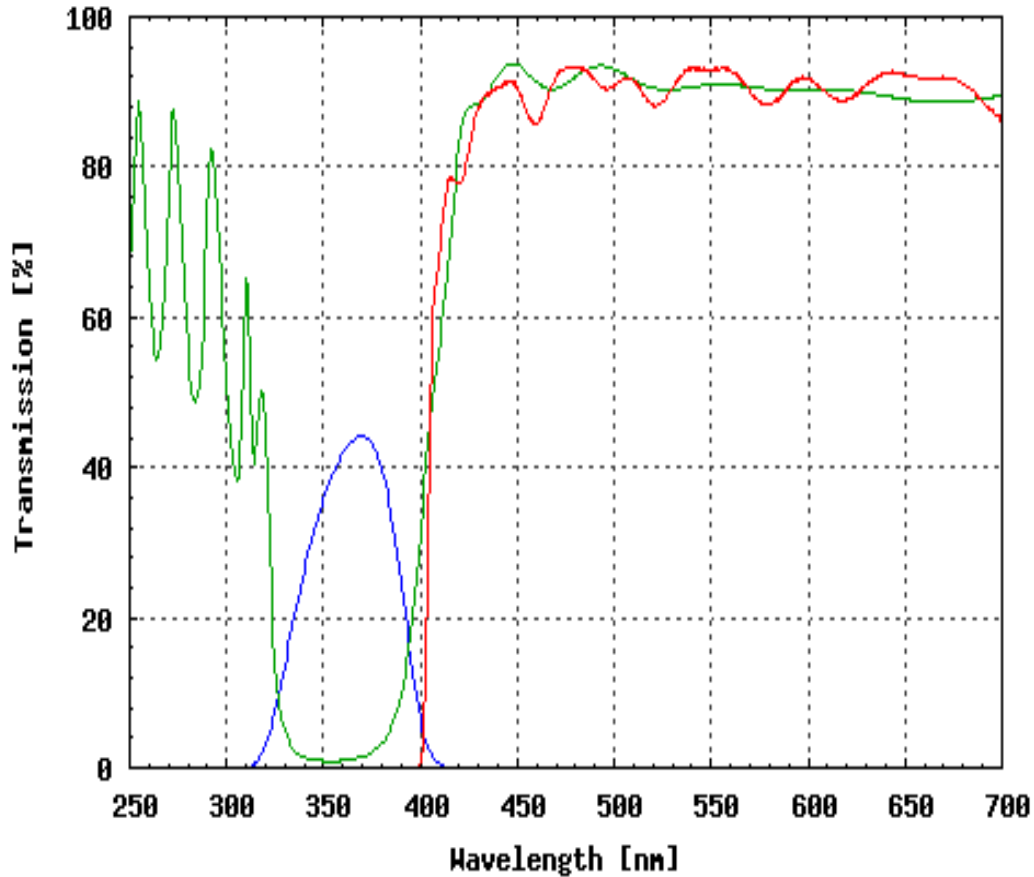

### FILTER SET COMPONENTS

**Exciter:** XF1005 (365WB50)

**Dichroic:** XF2001 (400DCLP)

**Emitter:** XF3097 (400ALP)

### OTHER INFORMATION

Longpass blue set with a 365nm wide-band exciter. Brighter image with lower signal-to-noise than the XF06. Best for Hg source, not recommended for Xe.

See XF02 for a better longpass set choice for Xe.

Can be used for all Quantum Dot Conjugates, for simultaneous multi-colour viewing with Hg excitation