


## Omega Optical Filter Set XF100-2

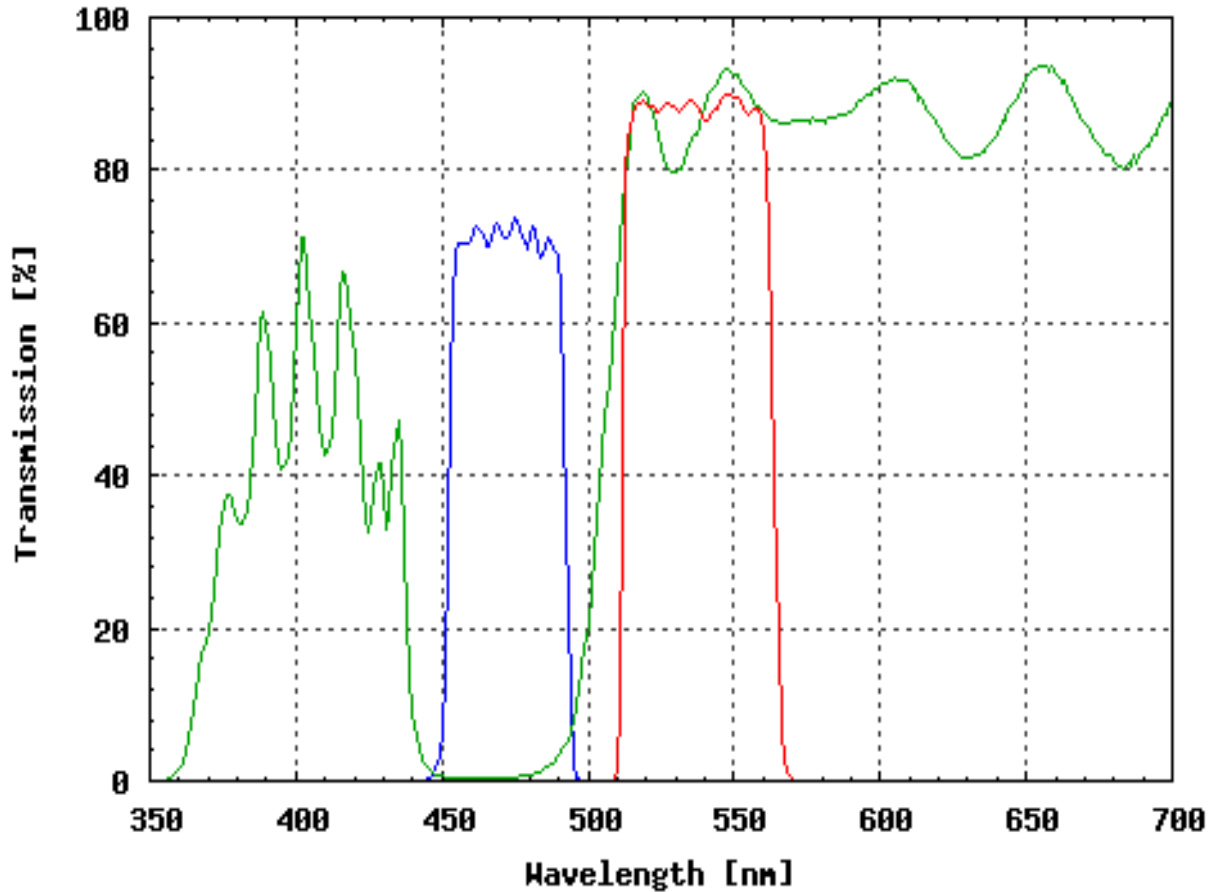

### FILTER SET COMPONENTS

Exciter:	XF1073 (475AF40)
Dichroic:	XF2010 (505DRLP)
Emitter:	XF3084 (535AF45)

### OTHER INFORMATION

This set XF100-2 option incorporates ALPHA technology. Highest quality green bandpass set for most green fluorophores. Higher signal-to-noise than the XF115-2 longpass equivalent and a better discriminating set than XF100-3. This set is outstanding as a GFP or FITC set, featuring a bright signal with a black background. The set is recommended for multiple labelling work, and will provide adequate spectral separation from dyes such as rhodamine, Cy3, or Cy5. See XF115 for a FITC or GFP longpass option.