


Omega Optical Filter Set XF101-2

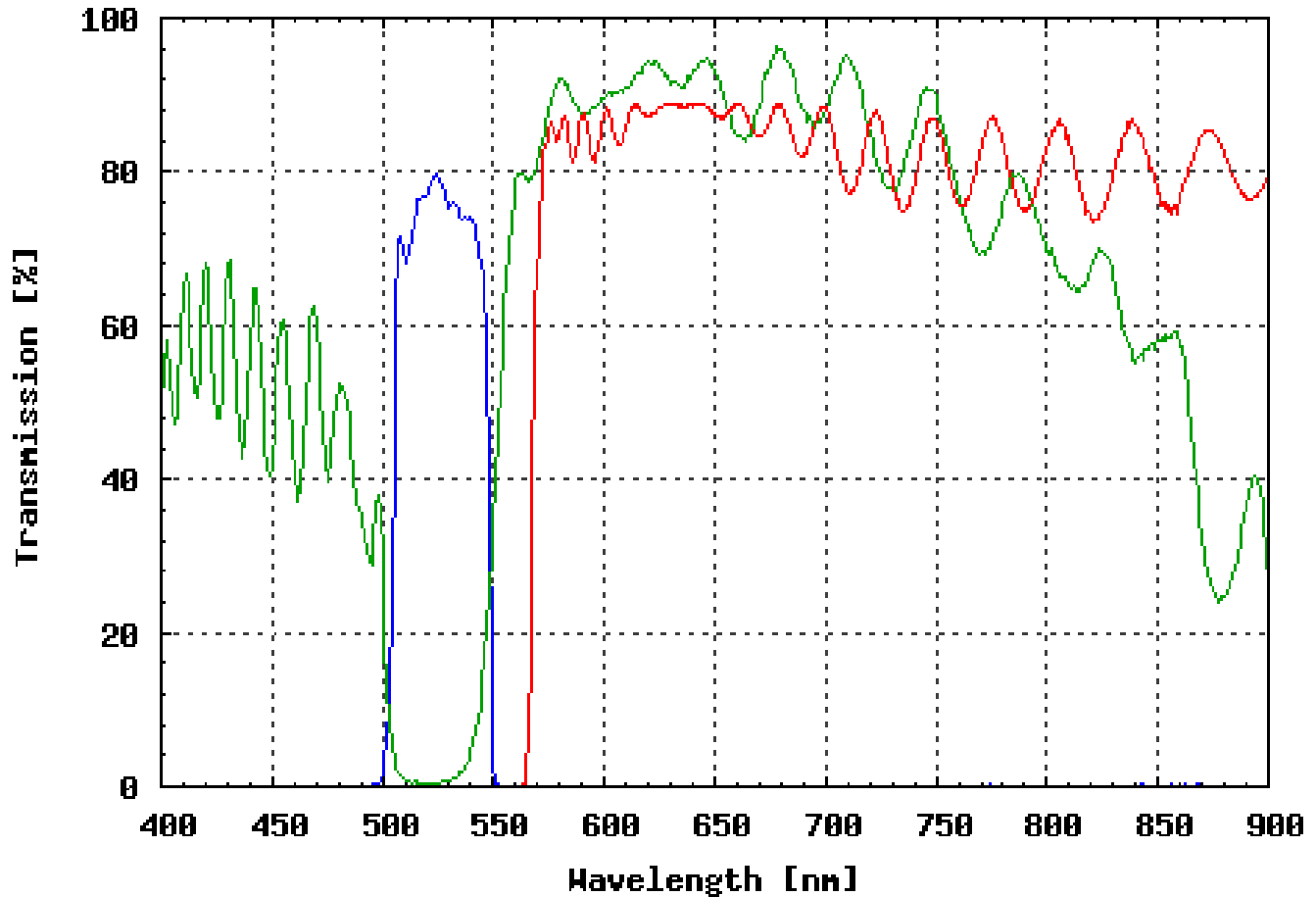

FILTER SET COMPONENTS

Exciter:	XF1074 (525AF45)
Dichroic:	XF2017 (560DRLP)
Emitter:	XF3085 (565ALP)

OTHER INFORMATION

Brighter emission with lower signal-to-noise than the XF108-2 bandpass equivalent. This set is excellent for working with rhodamine or Cy3, and with the longpass emission filter, this set is designed to collect maximal signal from all recommended dyes. However, we do not recommend this set if you intend to also use Cy5 in your experiments. See XF108 for a set suitable for that type of work.