


## Omega Optical Filter Set XF115-2

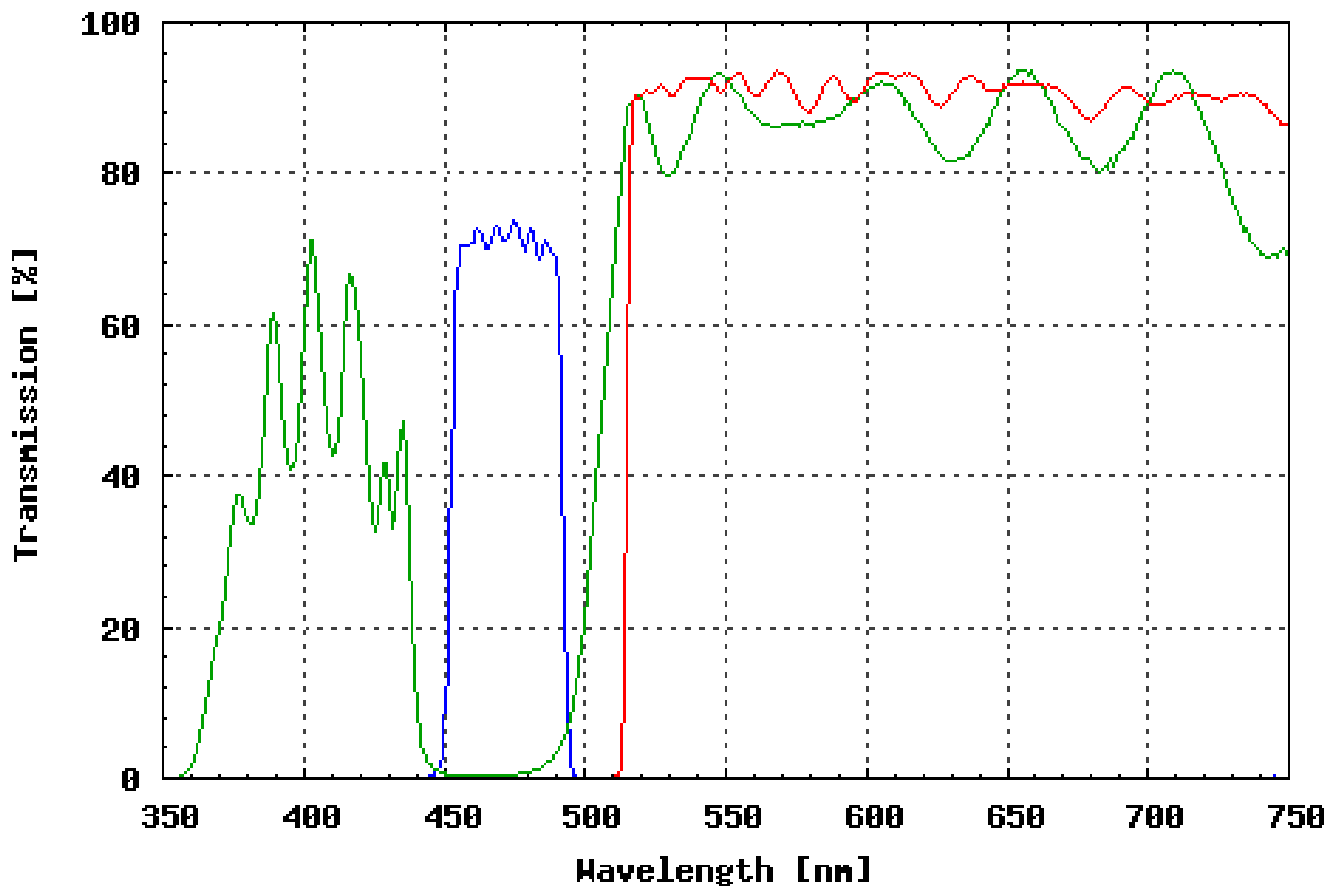

### FILTER SET COMPONENTS

<b>Exciter:</b>	XF1073 (475AF40)
<b>Dichroic:</b>	XF2010 (505DRLP)
<b>Emitter:</b>	XF3086 (510ALP)

### OTHER INFORMATION

Highest quality green longpass set for most green fluorophores. Brighter emission with lower signal-to-noise than the XF100-3bandpass equivalent. This set XF115-2 option incorporates ALPHA technology. The set has a better signal to noise ratio than our previous XF115 set. This set is outstanding as a GFP or FITC set, featuring a brighter signal than the XF100 - but also with a higher background. This set is recommended only if you are not going to work with dyes to the red, such as rhodamine, Cy3, or Cy5 - you will see bleedthrough from those dyes.