


Omega Optical Filter Set XF135-1

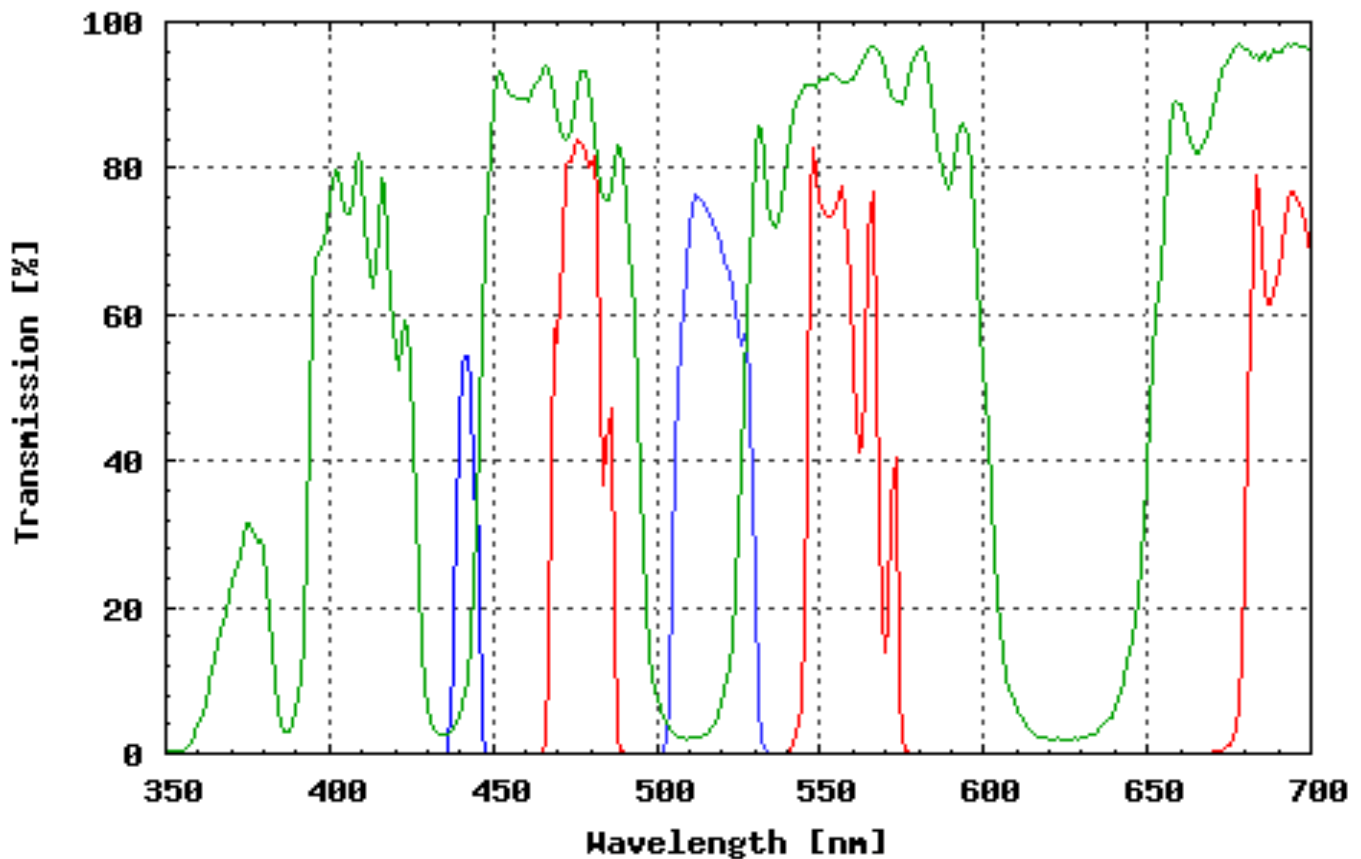

FILTER SET COMPONENTS

Exciter:	XF1079 (436DF10)
Exciter2:	XF1080 (510DF25)
Dichroic:	XF2065 (436-510DBDR)
Emitter:	XF3099 (475-550DBEM)

OTHER INFORMATION

A series of singleband excitors in a filter wheel or slider with a multi-band dichroic & emission filter housed in the microscope. This set is used to collect images for CFP-YFP with a PMT or a CCD camera. Letter designations are custom filter specific for the filter set being shipped & may not necessarily match the spectral characteristics exactly as displayed. Please note these excitors can not be ordered individually.