


Omega Optical Filter Set XF26

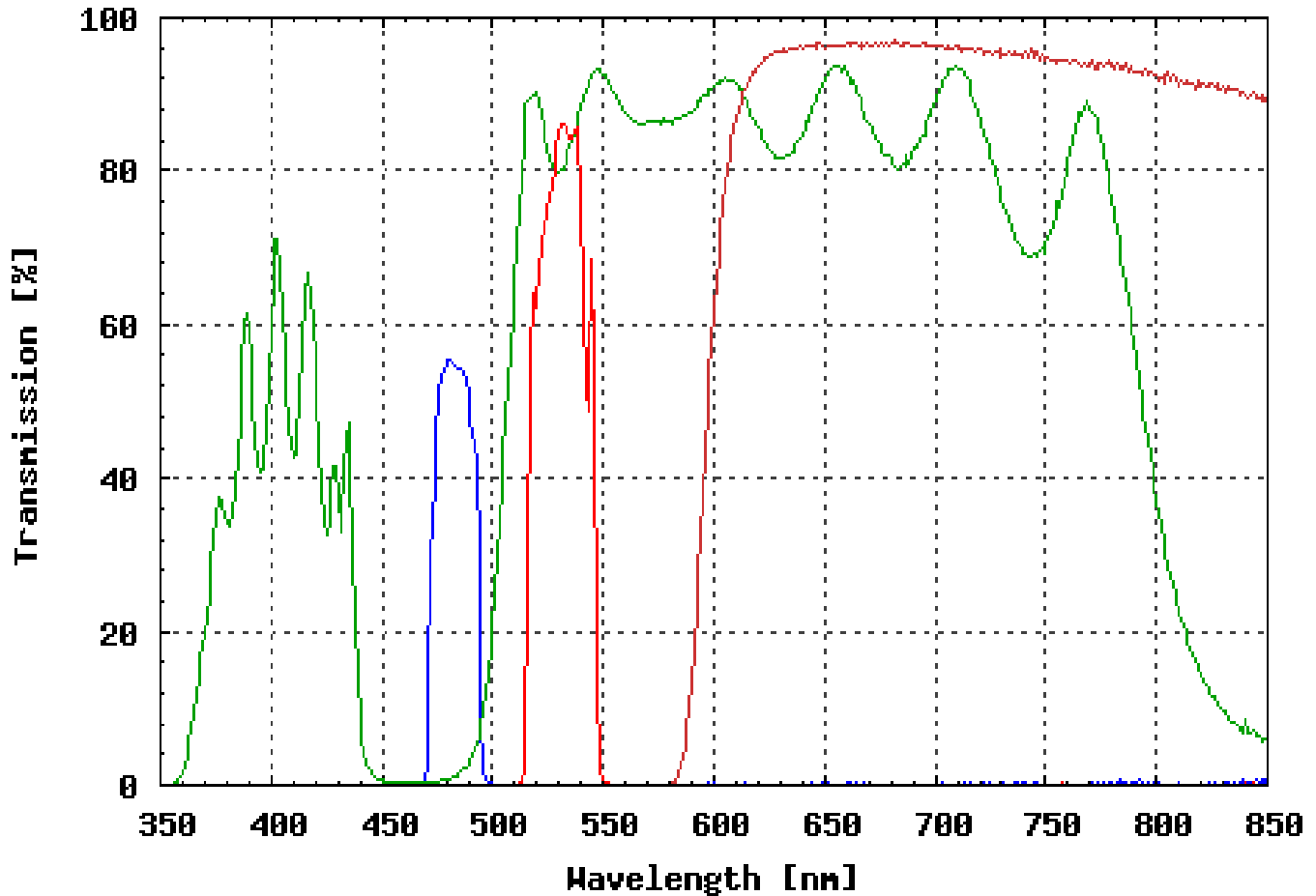

FILTER SET COMPONENTS

Exciter:	XF1015 (485DF22)
Dichroic:	XF2010 (505DRLP)
Emitter:	XF3016 (OG590)
Emitter2:	XF3017 (530DF30)

OTHER INFORMATION

This set is designed to allow the user to see one colour at a time (in contrast to the XF115, which allows simultaneous visualisation of both colours). Set is designed to work with an emission filter wheel. If an emission filter wheel is not available, one could configure two sets: Set 1 (XF1015, XF2010, XF3017) and Set 2 (XF1015, XF2010, XF3016). The user could then alternate between these sets to see the individual emission colours. Please contact us for further details or ordering information.

