


Omega Optical Filter Set XF50-1

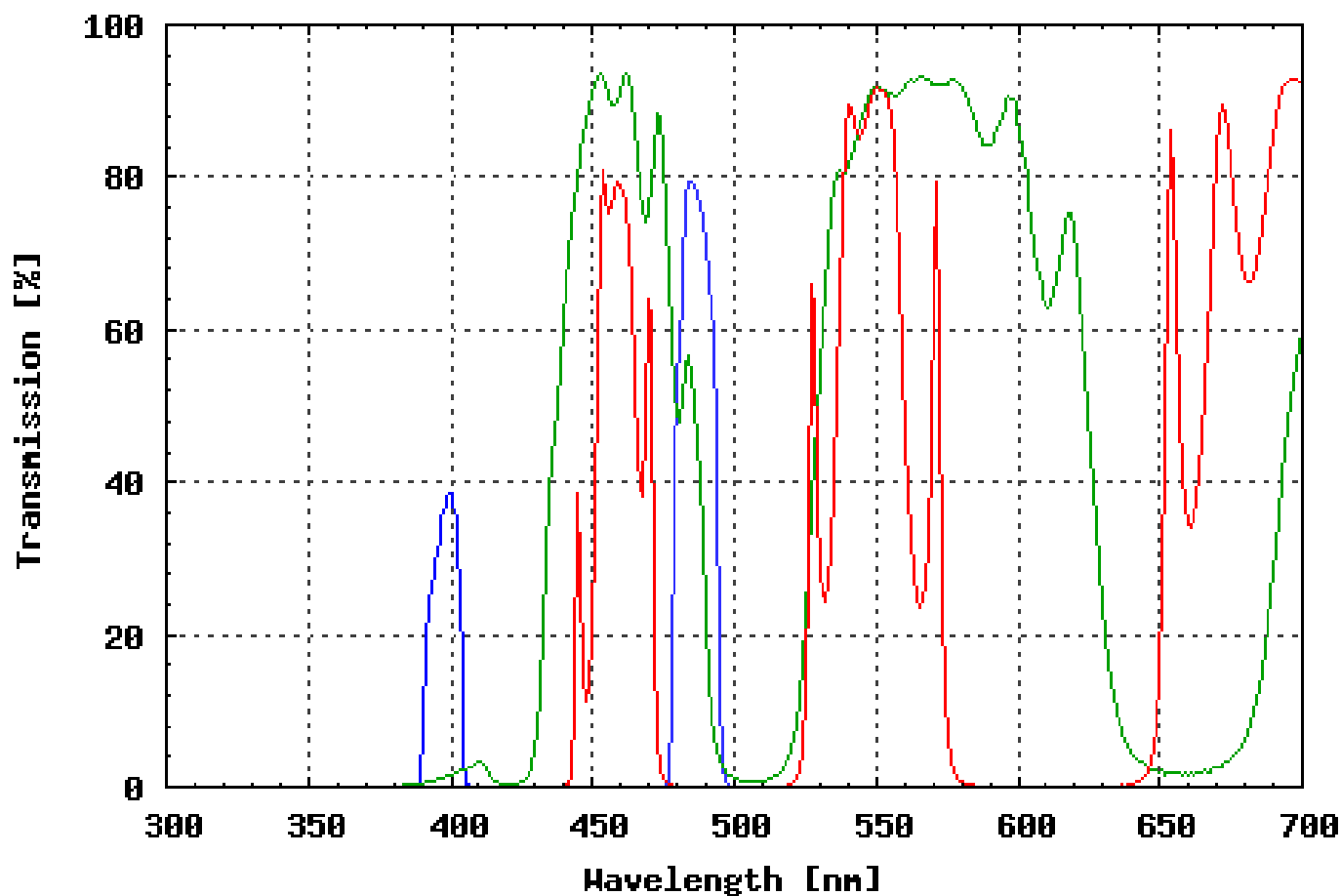

FILTER SET COMPONENTS

Exciter:	XF1006 (400DF15)
Exciter2:	XF1042 (485DF15)
Dichroic:	XF2041 (385-502DBDR)
Emitter:	XF3054 (460-550DBEM)

OTHER INFORMATION

A series of singleband exciters in a filter wheel or slider with a multi-bands dichroic and emission filter housed in the microscope. This set is used to collect images for DAPI-FITC or BFP-EGFP (enhanced) with a PMT or CCD camera. Letter designations are custom filter specific for the filter set being shipped and may not necessarily match the spectral characteristics exactly as displayed. Please note these exciters can not be ordered individually.