


Omega Optical Filter Set XF52-1

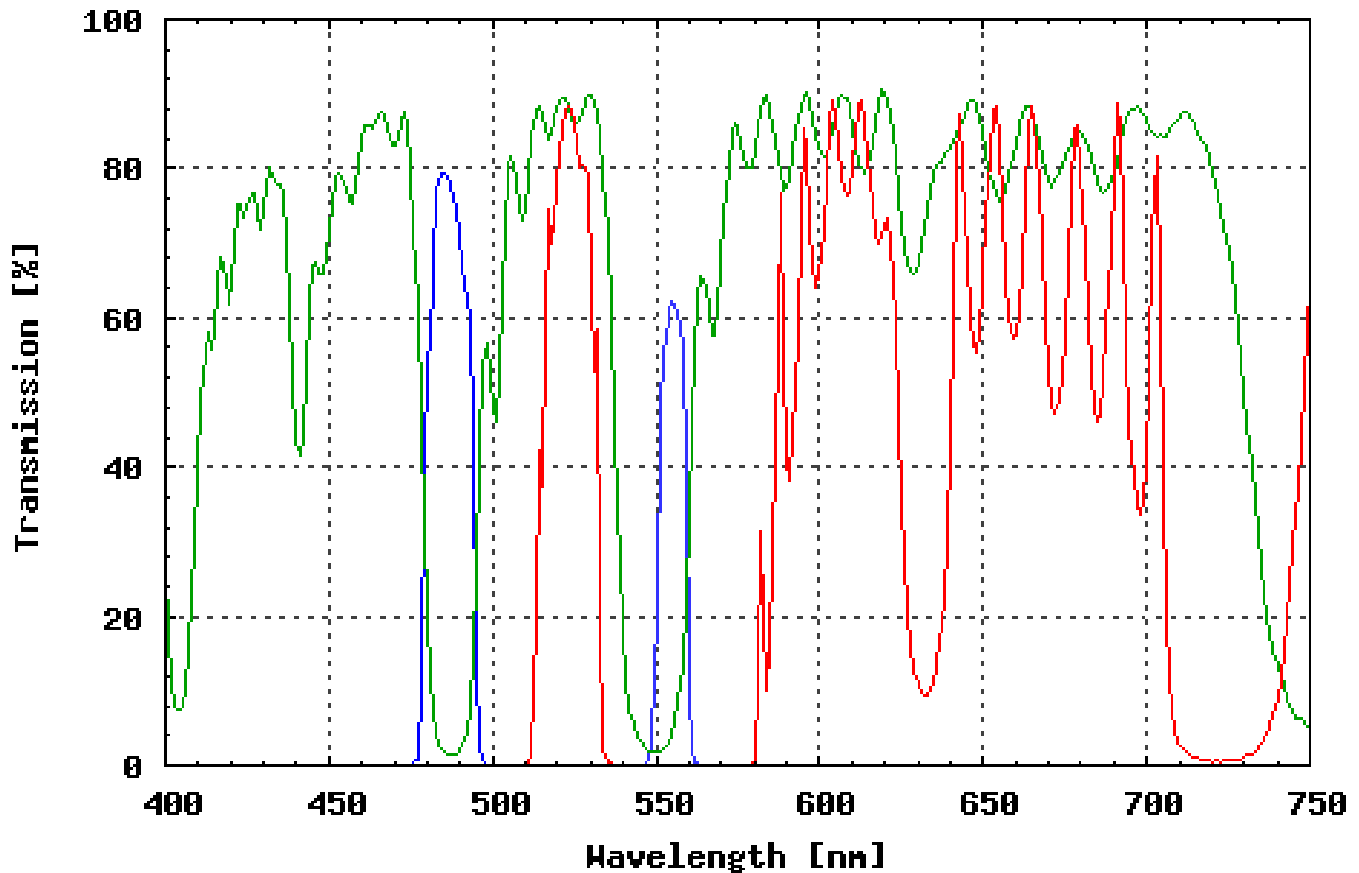

FILTER SET COMPONENTS

Exciter:	XF1042 (485DF15)
Exciter2:	XF1043 (555DF10)
Dichroic:	XF2043 (490-550DBDR)
Emitter:	XF3056 (520-580DBEM)

OTHER INFORMATION

A series of singleband exciters in a filter wheel or slider with a multi-band dichroic & emission filter housed in the microscope. This set is used to collect images for EGFP(enhanced)-DsRed2, FITC-TRITC or Cy2-Cy3 with a PMT or a CCD camera.. Letter designations are custom filter specific for the filter set being shipped & may not necessarily match the spectral characteristics exactly as displayed. Please note these exciters can not be ordered individually.