


Omega Optical Filter Set XF59-1

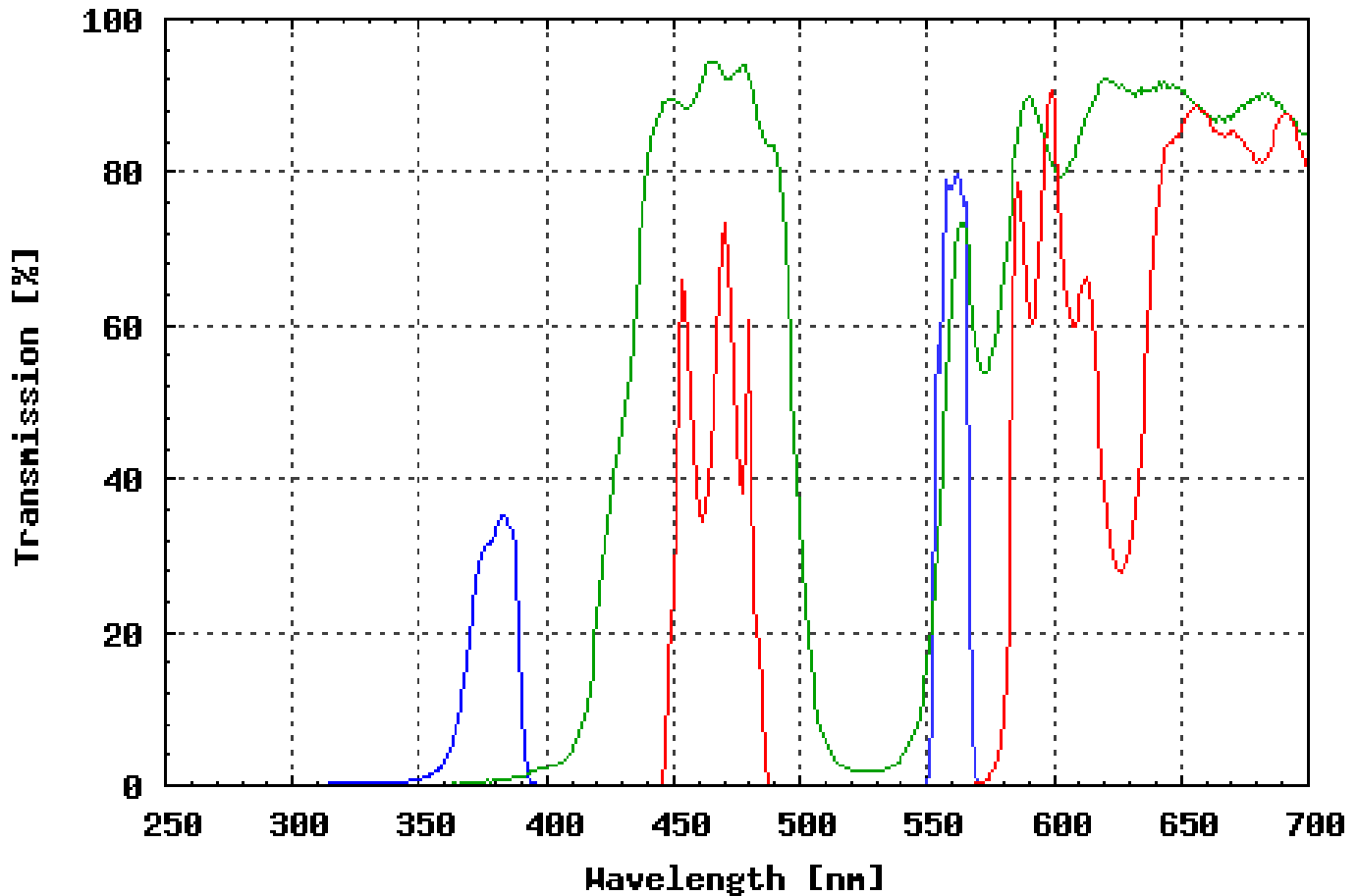

FILTER SET COMPONENTS

Exciter:	XF1094 (380AF15)
Exciter2:	XF1045 (560DF15)
Dichroic:	XF2047 (395-540DBDR)
Emission:	XF3060 (470-590DBEM)

OTHER INFORMATION

Pinkel Set for DAPI/TRITC

A series of singleband excitors in a filter wheel or slider with a multi-band dichroic & emission filter housed in the microscope. This set is used to collect images for DAPI-TRITC or BFP-DsRed2 with a PMT or a CCD camera..

Letter designations are custom filter specific for the filter set being shipped & may not necessarily match the spectral characteristics exactly as displayed.

Please note these excitors can not be ordered individually.