


## Omega Optical Filter Set XF68-1

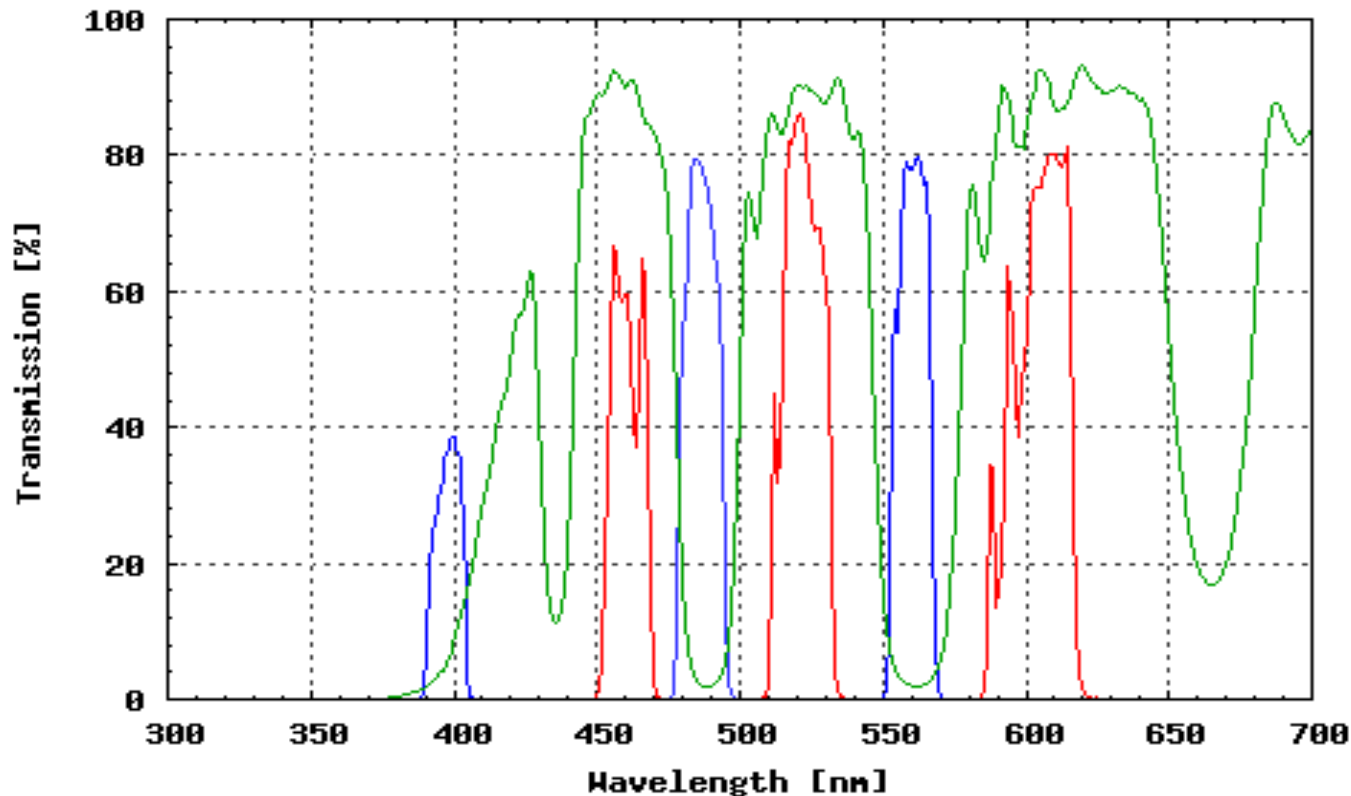

### FILTER SET COMPONENTS

|           |                          |
|-----------|--------------------------|
| Exciter:  | XF1006 (400DF15)         |
| Exciter2: | XF1042 (485DF15)         |
| Exciter3: | XF1045 (560DF15)         |
| Dichroic: | XF2050 (385-485-560TBDR) |
| Emitter:  | XF3063 (460-520-602TBEM) |

### OTHER INFORMATION

Pinkel Set for DAPI/FITC/TRITC. Could also be used for DAPI/Alexa Fluor 488/Alexa Fluor 546, DAPI/Cy2/Cy3, and other spectrally similar dyes. A series of singleband excitors in a filter wheel or slider with a multi-band dichroic and emission filter housed in the microscope. This set is used to collect images for DAPI-FITC-TRITC or DAPI-FITC-Cy3 with a PMT or a CCD camera. Letter designations are custom filter specific for the filter set being shipped & may not necessarily match the spectral characteristics exactly as displayed. Please note these excitors can not be ordered individually.