


Omega Optical Filter Set XF69-1

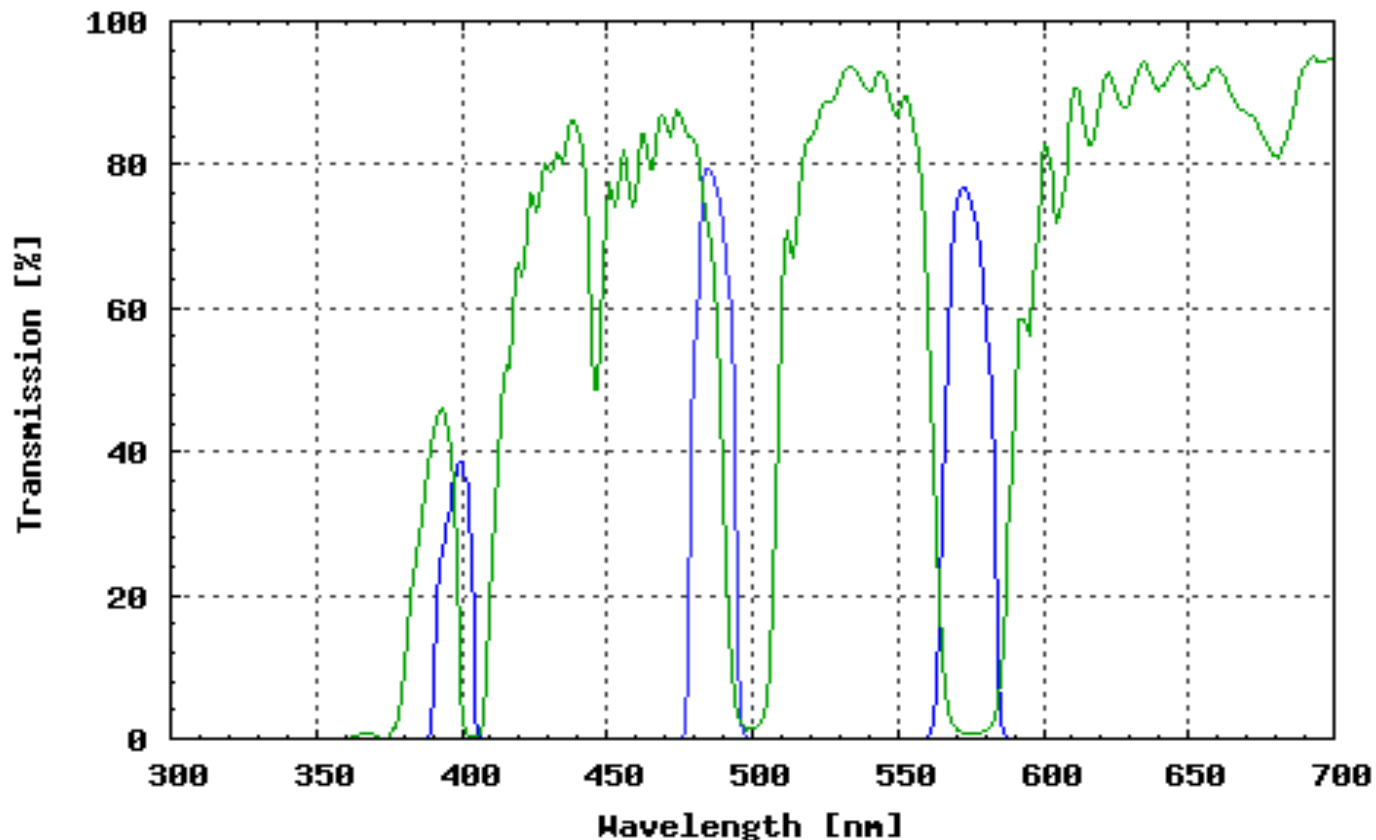

FILTER SET COMPONENTS

Exciter:	XF1006 (400DF15)
Exciter2:	XF1042 (485DF15)
Exciter3:	XF1044 (575DF25)
Dichroic:	XF2051 (400-495-575TBDR)
Emitter:	XF3116 (470-530-620TBEM)

OTHER INFORMATION

Pinkel Set for DAPI/FITC/Propidium Iodide. A series of singleband excitors in a filter wheel or slider with a multi-band dichroic & emission filter housed in the microscope. This set is used to collect images for DAPI-FITC-Propidium Iodide with a PMT or a CCD camera.

Letter designations are custom filter specific for the filter set being shipped & may not necessarily match the spectral characteristics exactly as displayed. Please note these excitors can not be ordered individually.