


Omega Optical Filter Set XF79-2

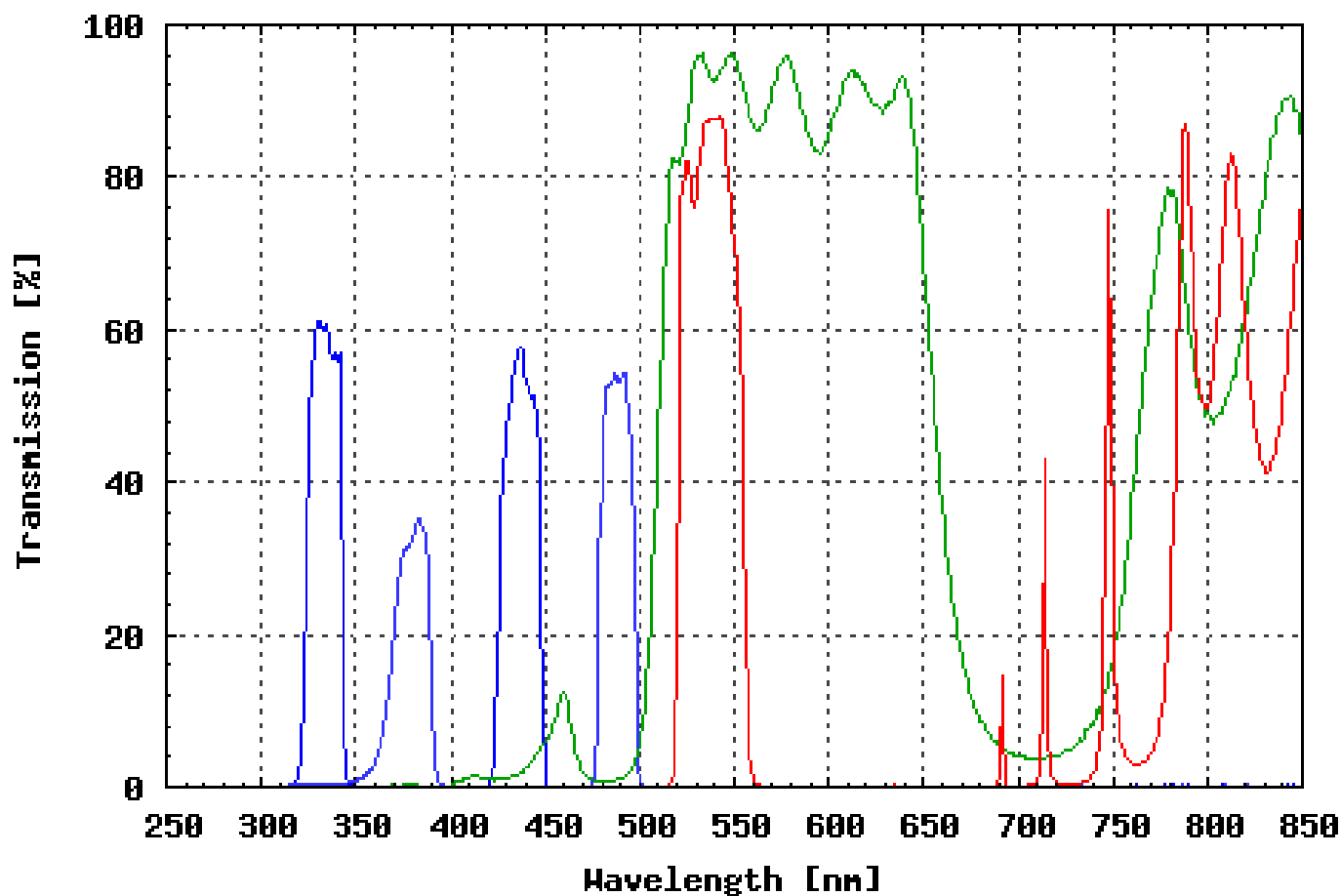

FILTER SET COMPONENTS

Exciter: XF1011 (490DF20)

Exciter: XF1071 (440AF21)

Exciter: XF1093 (340AF15)

Exciter: XF1094 (380AF15)

Dichroic: XF2031 (505DRLPXR)

Emitter: XF3007 (535DF35)

OTHER INFORMATION

Fura-2/BCECF Dual Ratio Imaging Set - This set is designed to allow for simultaneous ratio imaging of BCECF and Fura-2. If you are working with a monochromator, order the dichroic and emission filter only.

