


Omega Optical Filter Set XF93-1

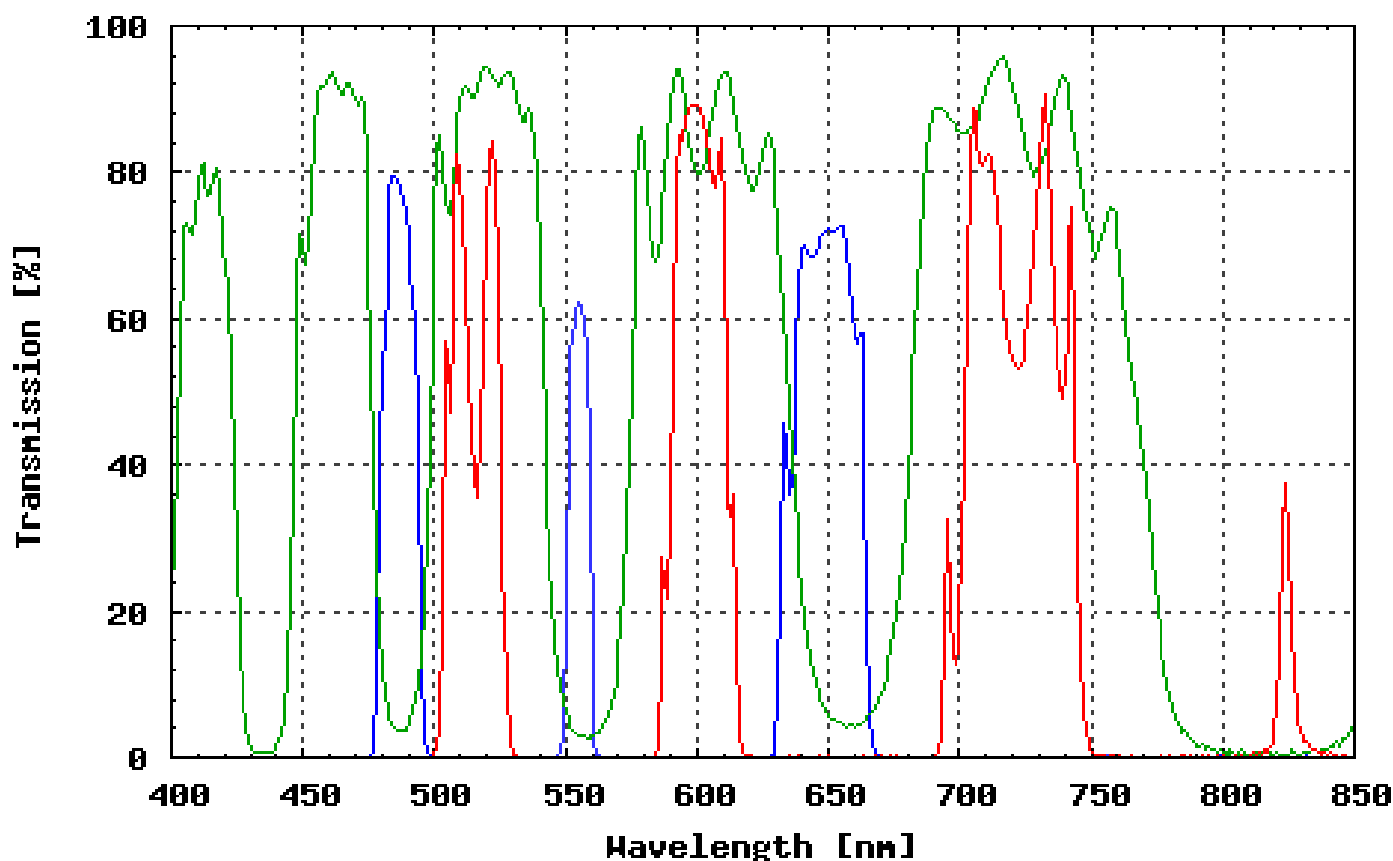

| <u>FILTER SET COMPONENTS</u> | |
|-------------------------------------|--------------------------|
| Exciter: | XF1042 (485DF15) |
| Exciter | XF1043 (555DF10) |
| Exciter | XF1046 (655DF30) |
| Dichroic: | XF2054 (485-555-650TBDR) |
| Emitter: | XF3067 (515-600-730TBEM) |

OTHER INFORMATION

A series of singleband exciters in a filter wheel or slider with a multi-band dichroic & emission filter housed in the microscope. This set is used to collect images for FITC-Cy3-Cy5 or FITC-TRITC-Cy5 with a PMT or a CCD camera. Letter designations are custom filter specific for the filter set being shipped & may not necessarily match the spectral characteristics exactly as displayed. Please note these exciters can not be ordered individually.